




1-800-339-3322 www.capemedical.com

PATIENT GUIDE IN USE OF COMPRESSOR/NEBULIZER

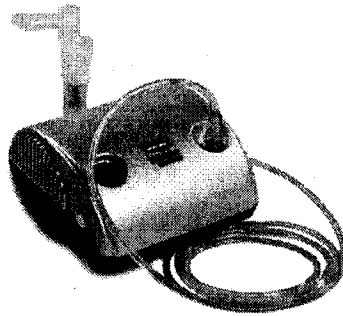

OPERATING INSTRUCTIONS

1. Before initial operation, nebulizer should be cleaned. Follow the recommended Cleaning Instructions outlined below.
2. **Do not run unit without filter.**
3. Place the device on a table so that the controls can be reached when you are seated.
4. Secure one end of the tubing to the compressor air outlet.
5. Connect the other end to the nebulizer.
6. Add liquid medication to the nebulizer according to your physician's prescription.
7. Plug into a 115V/60Hz wall outlet.
8. Insert the nebulizer into the tee.
9. Connect the mouthpiece to the tee.
10. Connect the reservoir tube to the tee.
11. Press the power switch ON to start the compressor.
12. Place mouthpiece between your teeth and close your mouth.
13. As aerosol begins to flow, **inhale deeply, slowly, and hold your breath for 10 seconds.** Then exhale slowly through the mouthpiece. Continue until medicine is gone. Treatment usually takes 10-15 minutes.

CLEANING INSTRUCTIONS

1. With proper maintenance your **reusable** nebulizer kit should last 6 months. Your **disposable** nebulizer kit should last 14 days. It is recommended that an extra kit be kept at all times.
2. After each treatment, rinse nebulizer and shake off excess water.
3. At the end of every day:
 - a. Remove the nebulizer, mouthpiece, and tee from the air tubing.
 - b. Disassemble nebulizer chamber.
 - c. Fill a plastic container or bowl with hot water and mild detergent solution.
 - d. Thoroughly clean the nebulizer, mouthpiece and tee in hot water and mild detergent solution.
 - e. Remove from the water/detergent solution and rinse with clear hot tap water.
 - f. Air dry thoroughly on a clean paper towel.
4. Once a week:
 - a. Fill a plastic container or bowl with hot water and vinegar solution (1 part vinegar to 3 parts water)
 - b. Soak in hot water/vinegar solution for 30 minutes.
 - c. Rinse with hot tap water again and air dry thoroughly.
 - d. Keep outer surface of tubing dust-free by wiping regularly with a damp cloth.

CLEANING OR REPLACING THE FOAM INTAKE FILTER

1. Clean the filter once a month or sooner if filter discolors.
2. Unplug unit.
3. Open cover of storage compartment to gain access to filter.
4. Remove filter and clean with soap and water. Let dry thoroughly before use.
5. Replace filter if it becomes clogged, torn or has a worn appearance.

SAFETY PRECAUTIONS

1. Do not remove outer case of this unit. All disassembly and maintenance of this unit **MUST** be done by a qualified service technician.
2. Unplug the unit before cleaning. **DO NOT** submerge in water to clean. Clean outer case with damp cloth.
3. Always unplug this product immediately after using.
4. Do not place or store unit where it can fall or be pulled into a tub or sink.
5. Close supervision is necessary when this product is used by, on, or near children and physically challenged individuals.
6. Never operate this product with a damaged cord or plug.
7. Keep the cord away from heated or hot surfaces.
8. Never drop or insert any object into any opening.
9. Never use while sleepy or drowsy.
10. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free from lint, hair and the like. Blocked openings may cause the unit to shut down.
11. This unit is oil-less. **DO NOT** lubricate.