




1-800-339-3322 – www.capemedical.com

HOMEFILL INSTRUCTIONS

1. Turn on oxygen concentrator (bottom machine).
2. Assure that concentrator is set no higher than 3 lpm.
3. Remove black plastic cap from fill connector on tank.
4. Set tank to "OFF".
5. Attach tank to Homefill Compressor.
6. Line up fill connector on tank with port on Compressor.
7. Push down on tank until it engages.
8. Pull up on silver collar to secure fit.
9. Turn on Homefill Compressor (top machine).
10. "Wait" light will illuminate for 1 to 3 minutes while compressor warms up.
11. "Filling" light will illuminate during tank filling
12. "Full" light will illuminate when tank is full.
13. Turn off Homefill Compressor.
14. Hold on to tank.
15. Push down on the silver collar. Tank will release.
16. Remove tank from compressor.

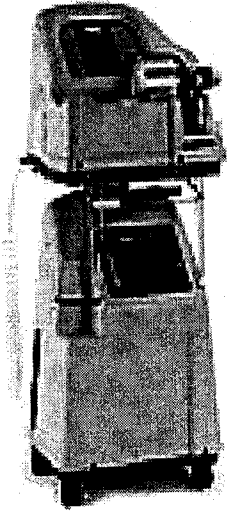

NORMAL CYLINDER TIME FOR EMPTY TANKS:

ML-4	60 minutes
ML-6	1 hour and 25 minutes
ML-9	2 hours and 20 minutes

